





JUMBOTM



Positioning chair **JUMBO™** is used for initial upright standing of children in the sitting position.

It is made of ecologic wood and therefore it is eco-friendly.

Despite its simplicity, it perfectly fulfils its function. It is safe for spastic as well as flaccid children with decreased muscle tone. Adjustable side supports help in correcting postural problems and, moreover, this chair may be beneficial regarding orthopaedic disorders.

It has a safety vest which is equipped with innovative "hold & pull" system thanks to which the child is both quick and precisely secured in this chair.

Each product is hand painted on water base paints. Size 3 thanks to the quick folding system of the footrest makes it easy to place and remove the person from the chair.








IT IS RECOMMENDED IN THE FOLLOWING CASES:

- Cerebral palsy various types and forms.
- Brain and spinal cord injuries with lower limbs paralysis and weakness, balance and motor coordination disorders with different etiology.
- Demyelinating disease (e.g. multiple sclerosis SM).
- Meningocele.
- Muscles diseases with weakness and paralysis (dystrophy, myasthenia gravis, myositis).
- Genetic disorders.
- Degeneration and neurological diseases.
- Other condition connected with muscles weakness, paralysis and musculoskeletal disorders.
- After head and spinal cord injuries and spine trauma with paralysis and paresis.

STANDARD EQUIPMENT:

Frame with castors and brakes | Foot rest with feet stabilizing belts | Seat | Seat cushion | Efficient Tilt-in-space adjustment
The use in the product of an innovative four-bar linkage, supported with the gas spring, greatly improves smooth Tilt-in-space adjustment. | Backrest

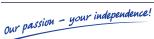



JUMBOTM



ADVANTAGES:



Stable structure

The device ensures adequate stability and safety of Patient.



Easy to use

Intuitive way of adjustment and easy assembly.



Eco-friendly materialS

We use only human and environment friendly materials.



Spacer upholstery

The device is equipped with breathable, hypoallergenic 3D spacer fabric. Now you can be sure that the skin stays dry for longer!





Hold and pull©

Our harnesses and pelvic belts are equipped with Hold and pull system that provides fast and comfortable support your child in the device.



Independently adjustable side supports ensure proper trunk stabilization.



Contoured abduction block ensures proper sitting position. Abduction block can be removed.



Due to available adjustments you are able to set the desired seat depth as well as tilt angle of the backrest.



PRODUCT IS AVAILABLE IN TWO OPTIONS:



Modern version



Classic version



JUMBOTM







Grass Size 1


Sea Size 2



Size 3



Optional colors available for special request. Production time ca. 14 days

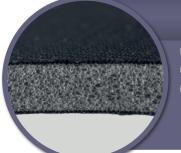

WATERPROOF UPHOLSTERY


- ☑ Lightweight
- ☑ Waterproof
- ☑ Easy to clean
- ☑ Ecological
- ☑ High durablity
- ☑ Easy decontamination



Green Grass Optional



Blue Sea Optional Optional colors available for special request. Production time ca. 14 days.

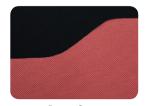

Brown Sugar Optional



Optional

Hot Red Optional



MEDICAL FABRIC

- ☑ Medical Fabric
- ☑ Breathable
- ☑ Waterproof
- ☑ Fire retardant*

easy to disinfection and because of thatthe product can be used by a many patients.Additional asset of the fabric is breathable and water-resistant. in hospitals.

*Dostępne For special request




SAN Blue Optional



SAN Skin Optional

Optional colors available for special request. Production time ca. 14 days.

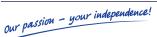


Akces-MED sp. z o.o. Jasionka 955B 36-002 JASIONKA

export@akces-med.com +48 (17) 864 04 77 www.akces-med.com







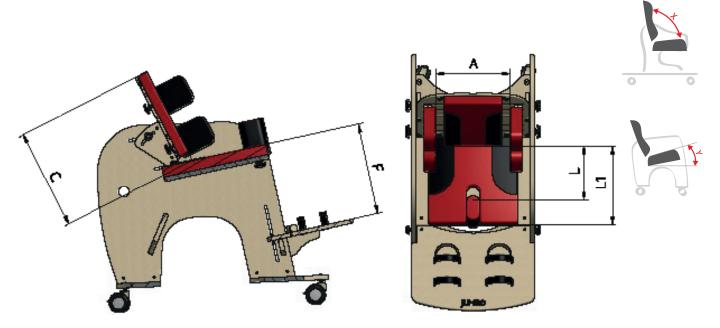

Symbol	Measurement	Unit	Size 1	Size 2	Size 3
А	Seat width	[cm]	35	39	46
L	Seat depth	[cm]	13-31	19-38	29-54
L1	Seat depth to the abduction block	[cm]	0-24	3-29	8-45
C	Backrest height	[cm]	34	38	55
Е	Backrest width	[cm]	17-27	22-32	22-35
F	Foot length adjustment	[cm]	18-35	24-45	28-45
	Footer lenght	[cm]	19	25	26
X	Seat tilt angle	[°]	90-110	90-110	90-110
Υ	Backrest tilt angle	[°]	10	10	10
	Max user weight	[kg]	30	45	75

Product dimensions and weight	Unit	Size 1	Size 2	Size 3
Width	[cm]	50	56	60
Lenght	[cm]	82	100	120
Height	[cm]	76	95	105
Weight	[kg]	14	20	24,5

Note: Measurement tolerance range is +/- 1 cm / $\frac{1}{2}$ in.



SŁONIKTM





Safety



Soft cover for tray

Soft cover protects the user from hitting the tray. Indicated for patients with impaired limbs coordination.



Patient lift belt

Robust and ergonomic design as well as high durability facilitates daily care of the patient.



Stabilization



Adjustable headrest

Headrest adjustable in three dimensions. Due to its special design, it can be precisely adjusted to the user's needs.



Extended headrest

For patients with increased head circumference. It provides better side stabilization.



Head supporting belt

Attachable to the headrest by use of Velcro straps. Used for people with limited control of the head to prevent its falling.



Thigh abduction belts

They can be easily attached to the chair by use of Velcro straps. Their use allows to keep legs in abduction and also prevents the patient from slipping off of the device.



H vest

It acts as a trunk stabilizer. HOLD&PULL© system makes it easier to adjust the vest to the patient body. You can quickly detach the vest if necessary.



6 point vest

Wide range adjustable vest that acts as a trunk stabilizer. HOLD&PULL© system makes it easier to adjust the vest to the patient body.



SŁONIKTM





Comfort



Tray for manual therapy

The tray allows active exercises assisting upper limbs with a particular focus on palm function.



Tray

Made of wood. It supports the therapy and allows you to play, eat, and work.



Cotton cover for headrest

It helps to keep the headrest clean. Easy to assemble/disassemble, it can be washed.



Seat cushion - Elastico

It precisely adapts to the patient's body. Made of memory foam covered with breathable and spatial fabric. It is indicated in anti-bedsore prophylaxis.



Backrest cushion - Elastico

It precisely adapts to the patient's body. Made of dense memory foam covered with breathable and spatial fabric. It is indicated in anti-bedsore prophylaxis.



Seat cushion - profiled thighs

Anatomical notch on the thighs improves sitting comfort.



Seat cushion - Abduction block

Shaped abduction block provides with correct seating position. Can be easily removed



Tango caster with a diameter of 100mm

Noiseless, easy to maintain, does not leave marks, does not corrode.



Tray for non verbal communication

It is a very important element in the rehabilitation of patients with speech dysfunction and, moreover, it makes it easier to perform daily activities. The set consists of 38 cards that may be put in random order into 24 places of the tray. It is covered with transparent and aesthetic plate "ade by Plexiglas®. Parents and child can create their own symbols that they believe are significant for their communication.



Storage



Oxygen tank case

Oxygen tank case can be attached to the stroller by use of Velcro straps with adjustable length. Made to measure to perfectly fit the patient's oxygen tank.